II Bilancia Digitale

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Package contents

Scale
Alkaline batteries 2 x 1.5V, type AA
Instruction manual

With its ergonomic and compact shape, the Chicco Digital Scale weighs baby from the first few days with a capacity ranging from 30 g to 20 kg. The 10 g resolution also allows the measurement of small variations in weight; a sophisticated anti-vibration program stabilises weighing to measure accurately and quickly the weight of the baby even when he is moving.

Legend (Fig. A)

1. Scale weighing plate
2. LCD display
3. On/off button 
4. Tare function key 
5. Battery Compartment

WARNINGS

- This scale is only to be used for domestic purposes and cannot be used for public use or for commercial transactions.
- Always lift the scale when moving it. Do not drag or scrape it as it may become damaged or its functioning may be altered.
- Always use the scale with the utmost care.
- Never leave baby on the scale unattended!
- Only rest the scale on flat, stable surfaces.
- Never place weights on the scale exceeding its maximum capacity (20kg).
- Where the scale is not to be used for long periods, remove the batteries from the battery compartment.
- Keep the scale and the components of the packaging out of the reach of children.
- After having removed the appliance from its packaging, ensure that the scale is intact and that there is no visible damage. If in doubt, do not use it and only contact the retailer.
- Do not immerse the scale in water or other liquid, do not put it under a jet of water or allow water to get into it, or use it in the proximity of sources or accumulations of water: the scale could be damaged or its functioning could be altered.
- Product features may be subject to change without notice.

HOW TO USE

1. PREPARATION

Remove the battery compartment cover on the bottom of the scale (Fig. 1), insert the batteries supplied, taking care to observe the correct polarity as indicated on the inside of the compartment and replace the cover correctly. Position the scale on a horizontal smooth and stable surface without dragging it or scraping it along the surface (Fig. 2). Never place weights on the scale exceeding their maximum capacity (20kg).

2. SWITCHING ON AND DISPLAYING OF LAST WEIGHING

Press the on/off key , the display will indicate the symbols  for about 1 second. Then the value of the last weighing carried out will appear for about 2 seconds followed by the indication  (Fig.3). It is possible, thanks to the memory storage of the last weighing, to calculate the difference with the weighing in progress: simply make a note of the weight that will appear before performing the next weighing.


2.2 WEIGHT MEASUREMENT

When the scale display indicates , it is possible to place a towel on the weighing plate and then press the TARE key  (Fig. 4). Check that the display indicates  then place baby in the centre of the scale weighing plate (Fig. 5).

ATTENTION: the scale is equipped with an automatic device to stabilise the weight to eliminate the effects of measuring interference caused by the inevitable movements of baby. The waiting time may vary depending on the movements of baby. When the weight has stabilised, the display will flash 3 times, and after the third flash the weight of the baby will appear (Fig. 6). The scale will automatically switch off after 30 seconds (Fig. 7).

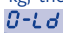ATTENTION: with weights greater than 20 kg, the scale will show the following error message  (Fig. 8).

3. FLAT BATTERIES

The display of the symbols  , indicates the flat status of the batteries which should be replaced.

3.1 REMOVAL, REPLACEMENT AND INSERTION OF THE BATTERIES

Remove the battery compartment cover on the bottom of the scale (Fig. 9) and remove the batteries. Insert the new batteries ensuring that they are of the same type or equivalent to the one indicated, being careful to observe correct polarity as indicated on the inside of the compartment; replace the cover correctly.

- The batteries must only be replaced by an adult.
- Use alkaline batteries equal or equivalent to those indicated in this manual.
- Do not mix different types of batteries or flat batteries with new ones.
- Do not leave batteries or any tools within the reach of children.
- Do not short-circuit the battery poles.
- In case of prolonged non-use of the product, always remove the batteries.
- Always remove flat batteries from the product to avoid any loss of fluid that can damage it.
- Dispose of flat batteries by recycling, do not throw away, do not dispose of in fire or in municipal waste.
- If the batteries are leaking liquid, replace them immediately, making sure to clean the battery compartment; thoroughly wash your hands if they have come into contact with the liquid.
- Do not attempt to recharge non-rechargeable batteries: they may explode.
- Using rechargeable batteries is not advisable as doing so may reduce the functionality of the product.
- The product is not designed to run on replaceable lithium batteries.

WARNING: improper use may create hazardous conditions.

4. CLEANING AND MAINTENANCE

Clean the scale by using only a cloth dampened with water, do not use abrasive products or chemicals. Do not soak or immerse the scale in water. When not in use, store the scale in a cool dry place protected from light and/or heat sources, avoiding placing objects on top of them.

TECHNICAL SPECIFICATIONS

Mod. 06746/A

REF. 00 005577 000 000

Size: mm 610 x 345 x 64 mm

Min. capacity 0.03 kg (30 g)

Max. capacity 20 kg (20.000 g)

Graduation: 10 g

Power supply: 2 x 1.5 V AA alkaline batteries

Use temperature range +10/ +40°C

Transport and storage temperature range -10+50°C

FUNCTIONS

Weight lock (Vibration dampers)

Automatic shut-off

Tare button 

Last weighing storage

Low batteries indicator

Overload indicator

LCD display

SYMBOLS LEGEND

 Device conforms to the essential requirements of applicable CE directives

 WARNING!

 Read the instructions manual.

 Instructions for use

DECLARATION OF CONFORMITY.

Artsana S.p.A. hereby declares that this Chicco brand scale model 06746/A complies with all the essential requirements and other provisions indicated in the EU Directives 2004/108/EC (electromagnetic compatibility) and 2011/65/EU (RoHS). Please feel free to contact Artsana S.p.A. Via S. Catelli 1 – 22070 Grandate (Co) Italy for a copy of this EC Declaration.

 **This product complies with the Directive 2012/19/EU.**
The crossed bin symbol on the appliance indicates that the product, at the end of its life, must be disposed of separately from domestic waste, either by taking it to a separate waste disposal site for electric and electronic appliances or by returning it to your dealer when you buy another similar appliance. The user is responsible for taking the appliance to a special waste disposal site at the end of its life. If the disused appliance is collected correctly as separate waste, it can be recycled, treated and disposed of ecologically; this avoids a negative impact on both the environment and health, and contributes towards the recycling of the product's materials. For further information regarding the waste disposal services available, contact your local waste disposal agency or the shop where you bought the appliance.

 **THIS PRODUCT COMPLIES WITH EU DIRECTIVE 2006/66/EC**

The crossed bin symbol on the batteries indicates that, at the end of their life, they must be disposed of separately from domestic waste, either by taking them to a separate waste disposal site for batteries or by returning them to your dealer when you buy similar rechargeable or non-rechargeable batteries. The chemical symbols Hg, Cd, Pb, printed under the crossed bin symbol, indicate the type of substance contained

in the batteries: Hg=Mercury, Cd=Cadmium, Pb=Lead. The user is responsible for taking the batteries to a special waste disposal site at the end of their life, so that they can be treated and recycled. If the spent batteries are collected correctly as separate waste, they can be recycled, treated and disposed of ecologically; this avoids a negative impact on both the environment and human health, and contributes towards the recycling of the batteries' substances. Non-compliance with the norms on battery disposal damages the environment and human health. For further information regarding the waste disposal services available, contact your local waste disposal agency or the shop where you bought the batteries.

ARTSANA reserves the right to make changes to this instruction manual at any time and without warning. All reproduction, transmission, transcription or translation into other languages of this manual, or part of the same, is strictly prohibited without prior written authorisation by ARTSANA.

Warranty: This product is guaranteed for 2 years from the purchasing date for manufacturing defects.

Artsana warrants for a period of 2 years, from the purchasing date, that it will solve any defects due to faulty material or workmanship free of charge by repairing or replacing the scale. If you decide to use the guarantee, please remember to also send the purchase receipt when returning the weighing scales.

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